

dx compact **comfort**



 **EUROFLEX**



dx compact^{comfort} - flexible, strong, comfortable and...compact

DX Compact with its small chassis, strong motors and large airfilled driving wheels ideal for those who are looking for a power wheelchair for both indoor- and outdoor use in one. It will easily take you straight from the kitchen / work desk out to a walk with the dog or a trip to the local store. The comfortable Comfort seat in leather that supports all your body makes you sit ergonomically all day in DX Compact. The seat unit is available in 3 different sizes with low or high back rest. The seat depth is adjustable to as much as 61 cm so even long legged can find a correct posture in DX Compact. Seat functions as back angle, seat tilt, seat lift and leg supports can be ordered power assisted. (power seat lift is standard) Thanks to its small size and the quick folding of the back you can easily transport the DX Compact in station wagon. Dx Compact is also crash tested and approved as a drivers seat in a car together with the carsafe lock from Autoadapt.



Front wheel driven with strong motors and big airfilled wheels gives good obstacle climbing ability. You can sit a long time in the Comfortsitsen without fatigue.




The small and slim chassis makes it easy to drive around at home and in narrow shops.




Quick folding of the back rest makes it easy to transport DX Compact in the back of car like a station wagon.



power seat lift

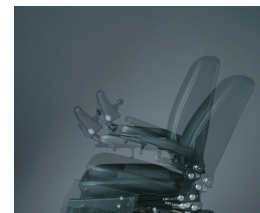

power tilt

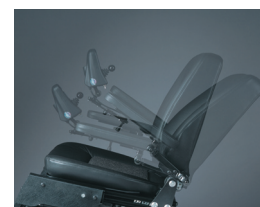

power back



power leg support

Standard features:

- chassis - dx compact - front wheel drive
- wheels - front airfilled - rear - solid
- Comfort seat unit- 45 cm, seat depth 40-61 cm
- seat lift - power
- back angle - manual - gas assisted
- seat tilt - manual - screw (onetime adjust))
- arm support - soft - fold back
- electronics - Dynamic - G90T

Accessories: (a selection)

- comfort seat unit - 40 cm, 53 cm - low back - 51 cm
- comfort seat unit - 45 cm - high back - 64 cm
- seat tilt - power
- back angle - power
- leg support - single foot plate
- leg support - divided manual or power
- head rest
- trunk support
- side support
- arm support - extra soft
- abductor
- table, tray

Technical facts:

Class	B, ISO122127
Obstacle climbing	65 mm
Drive	front wheels
Seat, sizes Comfort seat	40, 45, 53 x 51 cm, low back
Seat, sizes Comfort seat	45 x 64 cm, high back
Seat, depth -from back to seat front	40-61 cm
Seat, width- between arm support	46-56 cm
Seat, height- to seat top	45-63 cm,
Seat, tilt, power	-4° forward to +16° backward
Back, angle power	4° forward to -30° backward
Fabric upholstery	leather
Arm support, height	10-30 cm
Arm support, plate, size	31x8 cm
Lenght total - chassis	87 cm
Width total - chassis	62 cm
Weight total	120 kg
User weight max	125 kg
Turning area	117 cm
Driving time max	8,5 tim
Speed max	8 km/h
Electronics	Dynamic G90T

Eurovema AB reserves the right to make technical changes.



Eurovema
Creates mobility and ergonomics.

art.nr. BRO-900-ENG rev:2008-01