

miniflex ^{comfort}



miniflex^{comfort} - comfort in mini size

Miniflex Comfort is a power wheelchair for indoor use with minimal size and equipped with the Comfort seat, developed for users who spend a lot of time in their wheelchair. The popular and well known Miniflex can be ordered with this very comfortable ergonomical seat unit. Many customers appreciate the simple and reliable construction of the Miniflex, with it's flexibility and small size. The already good comfort is raised one more step with the Comfort seat. Miniflex is driving on the front wheels and the slightly larger driving wheels makes it easy to pass obstacles like threshold and the smallness makes it easy to get close to tables or working areas. The rear wheels are fixed in a joint to make good contact with the ground at all times.



In the Comfort seat you can sit a long time without fatigue.




The little and narrow chassis makes it easy to drive in your living area.



The quick release of the back support makes the wheel chair easy to transport in the trunk of a vehicle, e.g. station wagon.



power seatlift

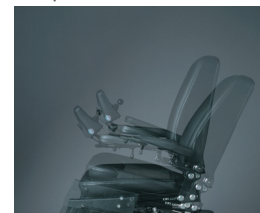

power tilt



power back



power leg support

Standard features:

- chassis - miniflex - front wheel drive
- wheels - frontwheel and rear wheel - filled
- comfort seat - 45 cm, seat depth 40-61 cm
- seat lift - power
- back recline - manual - gas assisted
- seat tilt - manual - screw (onetime adjust)
- arm support - soft - fold back
- elektronics - Dynamic - G90T

Accessories: (a selection)

- comfort seat unit - 40 cm, 53 cm -low back-51cm
- comfort seat unit - 45 cm - high back - 64 cm
- seat tilt - power
- back angle - power
- leg support - single foot plate
- leg support - divided manual or power
- head rest
- trunk support
- side support
- arm support - extra soft
- abductor
- table, tray

Technical data:

Class	A
Obstacles	max 30 mm
Driving	front wheels
Seat, width Comfortsits	40, 45, 53 cm
Back, height Comfortsits	51 low, 64 high cm
Seat, depth -from back to seat front	40-61 cm, adjustable
Seat, width- between arm support	46-56 cm, adjustable
Seat, height - to seat top	45-63 cm, alt. 50-75 cm
Sits, tilt, power	-4° forward to +16° backward
Back, recline	4° forward to -30° backward
Fabric	leather
Armrest, height	10-30 cm
Armrest, size	31x8 cm
Total lenght - incl. footrest	90 cm
Total lenght - chassis	73 cm
Total width - chassis	59 cm
Total weight	85 kg
Max user weight	125 kg
Turning space	120 cm
Max driving time	8,5 tim
Max speed	4 km/h
Electronics	Dynamic G90T

Eurovema AB reserves the right to make technical changes.



Eurovema
Creates mobility and ergonomics.

art.nr. BRO-500-03-ENG rev:2008-01