

Eurovema

Creates mobility and ergonomics.



VOLARIS
RA

Volaris RA

Volaris RA has been introduced for users who need extra support or have difficulty using a normal rollator. RA uses the same Comfort Frame™ used on all Volaris models. The rollator rolls smoothly over uneven ground. It is tested and approved for both indoor and outdoor use.

The rollator has arm supports that are adjustable in all directions, which provides safe and secure use. The soft wheels minimise vibrations in the handle, making it especially suitable for people with pain in their hands or those suffering from rheumatism. The cross-folding design means that it can stand stably on all four wheels when folded. It is easy to fold up and store in confined spaces.

The rollator is available in two frame heights and can be equipped with a range of accessories to meet user requirements. The Volaris RA is manufactured in Sweden from recyclable materials.

The rollator and accessories are tested and approved according to the applicable standards for walking aids,

SS-EN ISO 11199-3:2005.

Versions: standard and low.

Colours: natural or dark grey anodised aluminium.



Folded



Safety catch for height adjustment

Technical specification	RA	Low
Weight:	12 kg	11.9 kg
User weight:	125 kg	125 kg
Width:	64 cm	64 cm
Length:	80 cm	80 cm
Seat width:	47 cm	47 cm
Height:	100-125 cm	88-114 cm
Seat height:	46-57 cm	46-57 cm
Width when folded:	29 cm	29 cm
Wheel dimensions: (diameter x width)	20 x 3.2 cm	20 x 3.2 cm

Accessories



2490600 Tray



2497200 Drip stand



2431200 Bottle holder



2410230 Basket



2432200 Walking stick holder



2491000 Lamp



2494200 Soft seat



2461007 Single-handed brake




2460200 Parking brake



2441200 Oxygen cylinder holder



2413230 Textile basket holder

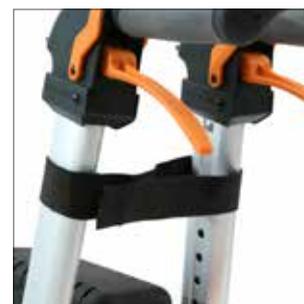

2495200 Transportation strap



2412200 Bag

Part no.	Designation
2490600	Tray Smart, grey
2497200	Drip stand
2431200	Bottle holder
2410230	Basket
2432200	Walking stick holder
2491000	Lamp
2494200	Soft seat
2461007	Single-handed brake RA

Part no.	Designation
2460200	Parking brake
2441200	Oxygen cylinder holder
2413230	Textile basket
2495200	Transportation strap
2412200	Bag
2412210-1	Bag*
2412220-1	Bag*

* Pattern and colour may vary



Eurovema

Creates mobility and ergonomics.

Eurovema Mobility AB

Baldersvägen 38, SE-332 35 Gislaved, Sweden

Tel. +46 (0) 371 390 100

Fax. +46 (0) 371 189 82

www.eurovema.se

info@eurovema.se

