

Light and flexible

Compact



Volaris S7 Compact the latest from Volaris family. Light and flexible walker for indoor with functions for everybody. Compact is also approved for outdoor use. With its 50 mm narrower framework makes it extra easy and maneuver in tight spaces and easier to get around. With frame construction

The Comfort Frame™ allows the rollator to remain stable on uneven surfaces. The frame flexes enabling the wheels to remain on the ground when going over bumps. The walker folds easily and repository standing with his x-folding design.




The seat and the handle can be easily adjusted in height allows it to accommodate multiple user groups. The walker is equipped with Smart handle. This means that there are no hanging wires that can hook into the subject, but everything is encapsulated in the frame. Both the process - and the parking brake is ergonomic designed and easily accessible and contribute to high safety.



Compact is made in Sweden of recyclable materials. Design and construction is patent pending. The walker is designed for existing standards for walking aids, ISO 11199-2:2005

Large accessory program to best meet the user's needs.

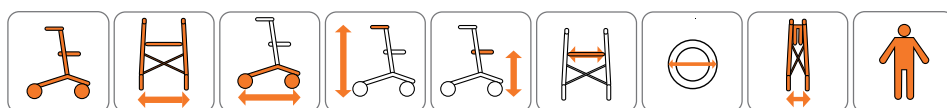

Compact	6,8 kg	53 cm	62 cm	73-88 cm	43-57 cm	36 cm	150x32 mm	27 cm	125 kg
Compact Low	6,8 kg	53 cm	62 cm	63-79 cm	43-47 cm	36 cm	150x32 mm	27 cm	125 kg

